

California and Oregon's Complete Residual Pesticide Analysis using a Shimadzu LCMS-8060

Jennifer Davis¹, Jeffrey H. Dahl¹, Christopher Gilles¹, Bob Clifford¹

¹ Shimadzu Scientific Instruments, Columbia, MD

1. Introduction

With the increase in medicinal and recreational cannabis legislation throughout the United States there is an emerging demand for pesticide testing on cannabis products. Currently each state is setting individual regulatory guidelines. This results in variation between the number of analytes tested and their required action levels; currently California regulates a total of 66 pesticides and Oregon regulates a total of 59 pesticides. To analyze these complete lists, laboratories commonly use both LCMS and GCMS because some compounds do not ionize well by ESI-LCMS. This study evaluates using both ESI and APCI ionization techniques to quantitate the complete California list using only LCMS. The resulting APCI-LCMS and ESI-LCMS MRM methods was tested in cannabis flower extract on a Shimadzu LCMS-8060. The LOQ determined for each pesticide was below the regulatory action level (Table 1).

2. Methods

A Shimadzu LCMS-8060 triple quadrupole mass spectrometer coupled with a Shimadzu Nexera X2 UHPLC system was employed for this evaluation. A total of 10 pesticides were analyzed by atmospheric pressure chemical ionization liquid chromatography mass spectrometry (APCI-LCMS). All other pesticides were analyzed by electro spray ionization liquid chromatography mass spectrometry (ESI-LCMS). Final method conditions can be seen in Figure 1 and 2. For each pesticide one to five MRM transitions were acquired. Separation was accomplished and retention times determined on column using neat standards prior to in-matrix evaluation. Matrix-matched calibration curves were prepared by serial dilution of spiked flower extract with blank flower extract and evaluated for each pesticide. The calibration set included multiple different concentrations, ranging from 1.0 ng/g to 2000 ng/g.

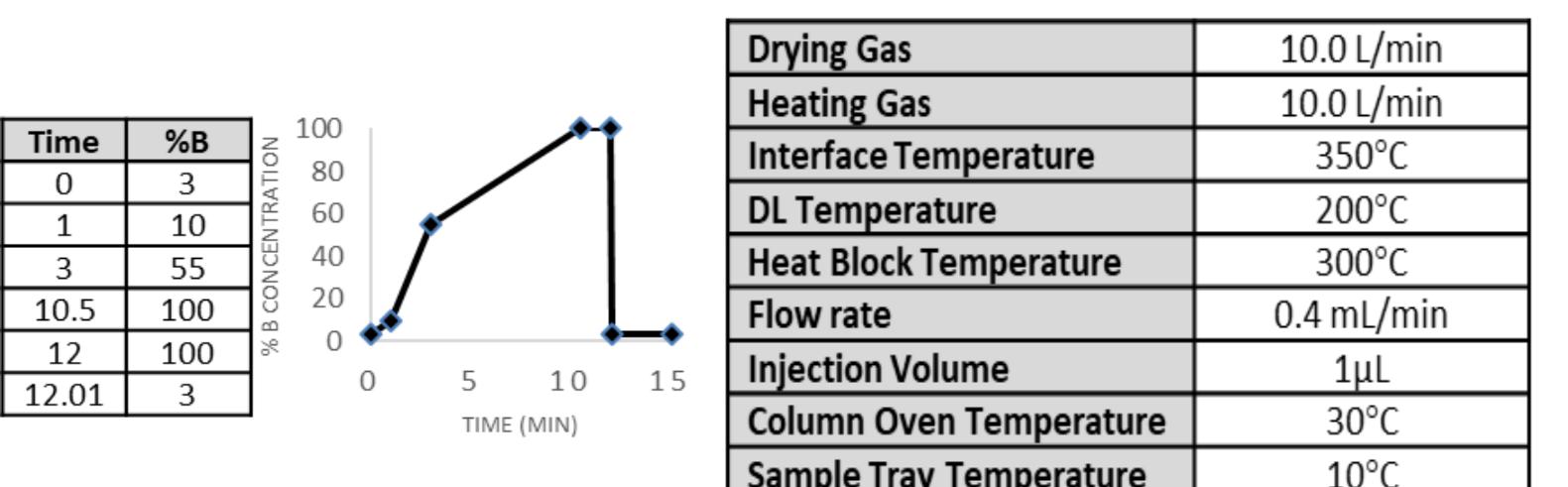


Figure 1 LCMS-8060 APCI Method Conditions

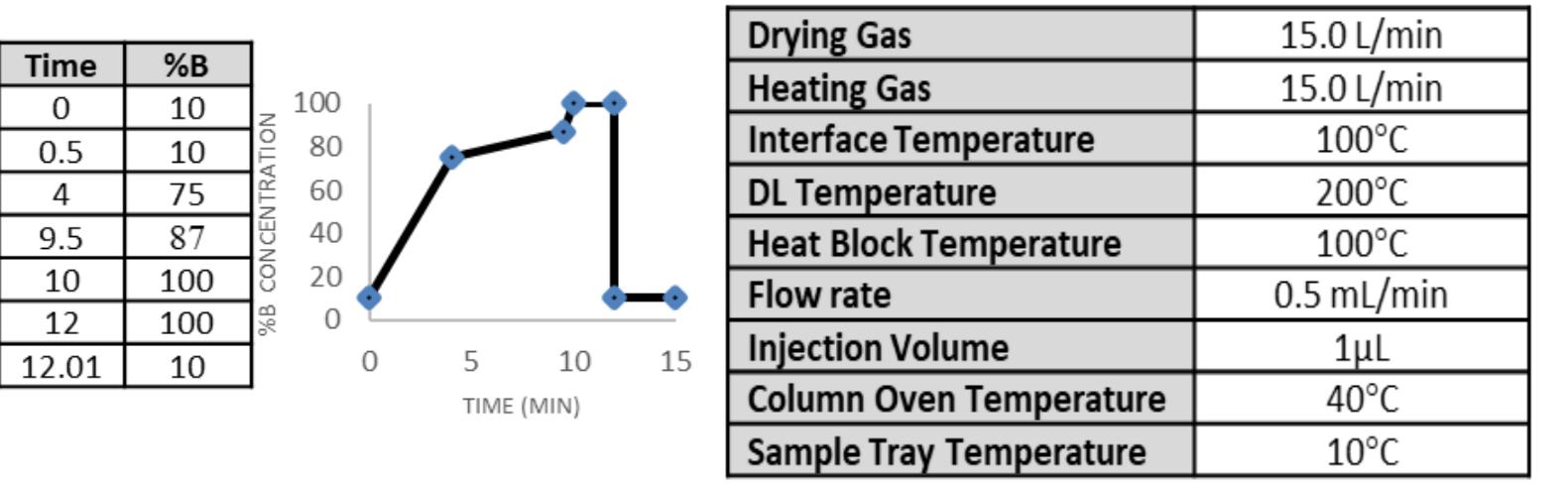


Figure 2 LCMS-8060 ESI Method Conditions

Table 1 LOQ determined for each Pesticide and the mode of Ionization

Residual Pesticide	LOQ (ng/g)	Method of Ionization	Residual Pesticide	LOQ (ng/g)	Method of Ionization
Abamectin (-)	31	ESI or APCI	Fluopyram (+)	2	ESI
Acephate (+)	20	ESI	Hexythiazox (+)	15	ESI
Acequonocyl (+/-)	15	ESI or APCI	Imazalil (+)	10	ESI
Acetamiprid (+)	<2	ESI	Imidacloprid (+)	4	ESI
Aldicarb (+)	<2	ESI	Kresoxim-methyl (+)	4	ESI
Allithrin (+)	50	ESI	Malathion (+)	2	ESI
Azoxystrobin (+)	4	ESI	Metalaxyl (+)	2	ESI
Bifenazate (+)	2	ESI	Methiocarb (+)	4	ESI
Bifenthrin (+)	4	ESI	Methomyl (+)	<2	ESI
Boscalid (+)	4	ESI	Methoprene (+)	50	ESI
Buprofezin (+)	<2	ESI	Methyl-parathion (-)	8	APCI
Captan (-)	8	APCI	Mevinphos (+)	4	ESI
Carbaryl (+)	10	ESI	MGK 264 (+)	15	APCI
Carbofuran (+)	<2	ESI	Myclobutanil (+)	10	ESI
Chlorantraniliprole (+)	2	ESI	Naled (+)	2	ESI
Chlordane (-)	15	APCI	Novaluron (-)	15	ESI
Chlorfenapyr (-)	8	APCI	Oxamyl (+)	2	ESI
Chlorpyrifos (+)	10	ESI	Paclobutrazol (+)	2	ESI
Clofentazine (+)	4	ESI	Pentachloronitrobenzene (PCNB) (-)	63	APCI
Clothianidin (+)	4	ESI	Permethrin (+)	10	ESI
Coumaphos (+)	4	ESI	Phenothrin (+)	10	ESI
Cyantraniliprole (-)	2	ESI	Phosmet (+)	10	ESI
Cyfluthrin (-)	15	APCI	Piperonyl butoxide (+)	5	ESI
Cyfluthrin (+)	500	ESI	Pirimicarb (+)	2	ESI
Cypermethrin (+)	60	ESI	Prallethrin (+)	10	ESI
Cyprodinil (+)	10	ESI	Propiconazole (+)	60	ESI
Daminozide (+)	15	ESI	Propoxur (+)	2	ESI
Deltamethrin (+)	30	ESI	Pyrethrosterin (+)	10	ESI
Diazinon (+)	<2	ESI	Pyrethrins (+)	100	ESI
Dichlorvos (DDVP) (+)	15	ESI or APCI	Pyridaben (+)	2	ESI
Dimethoate (+)	<2	ESI	Resmethrin (+)	35	ESI
Dimethomorph (+)	5	ESI	Spinetoram (+)	2	ESI
Dinotefuran (+)	2	ESI	Spinosad (+)	<2	ESI
Dodemorph (+)	4	ESI	Spirodiclofen (+)	10	ESI
Endosulfan-sulfate (-)	4	ESI	Spiromesifen (+)	20	ESI
Ethoprophos (+)	2	ESI	Spirotetramat (+)	2	ESI
Etofenprox (+)	4	ESI	Spiroxamine (+)	2	ESI
Etoxazole (+)	<2	ESI	Tebuconazole (+)	2	ESI
Fenheximid (+)	20	ESI	Tebufenozide (+)	5	ESI
Fenoxy carb (+)	2	ESI	Teflubenzuron (-)	15	ESI
Fenpyroximate (+)	10	ESI	Tetrachlorvinphos (+)	4	ESI
Fensulfotion (+)	5	ESI	Tetramethrin (+)	4	ESI
Fenthion (+)	100	ESI	Thiacloprid (+)	<2	ESI
Fenvalerate (+)	100	ESI	Thiamethoxam (+)	<2	ESI
Fipronil (-)	2	ESI	Thiophanate-methyl (+)	5	ESI
Flonicamid (-)	25	ESI	Trifloxystrobin (+)	<2	ESI
Fludioxonil (-)	2	ESI			

3. Results

Using the same sample, the entire California and Oregon residual pesticide list was analyzed. The APCI-LCMS method demonstrated accurate and precise trace-level quantitation in cannabis flower for 10 pesticides that are traditionally analyzed by GCMS and the ESI-LCMS method demonstrated the same robust and reproducibility for the remaining 84 pesticides. Both gradient methods were successfully used for chromatographic separation and identification of all 94 pesticides (Figure 3 and 4). The LOQ for each pesticide was below the California and Oregon action levels in cannabis, and precision and accuracy results were excellent. LOQs were determined for each pesticide using their corresponding retention time and a S/N calculation above 10:1.

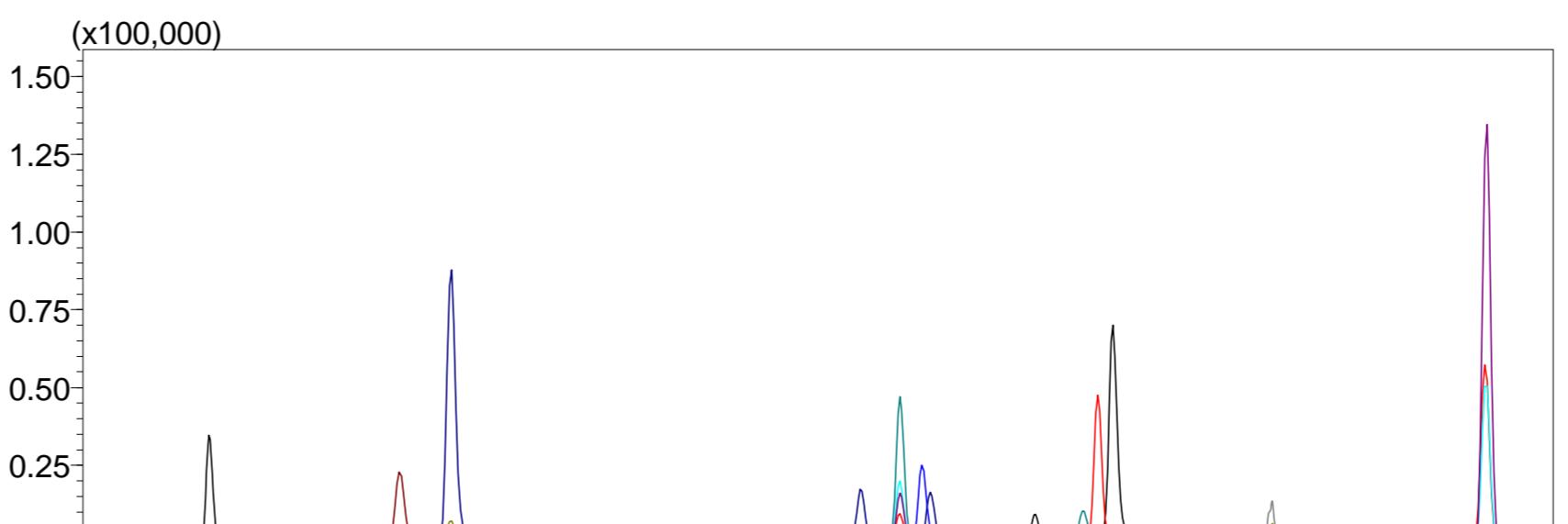


Figure 3 Representative Chromatogram for 10 pesticides using APCI-LCMS

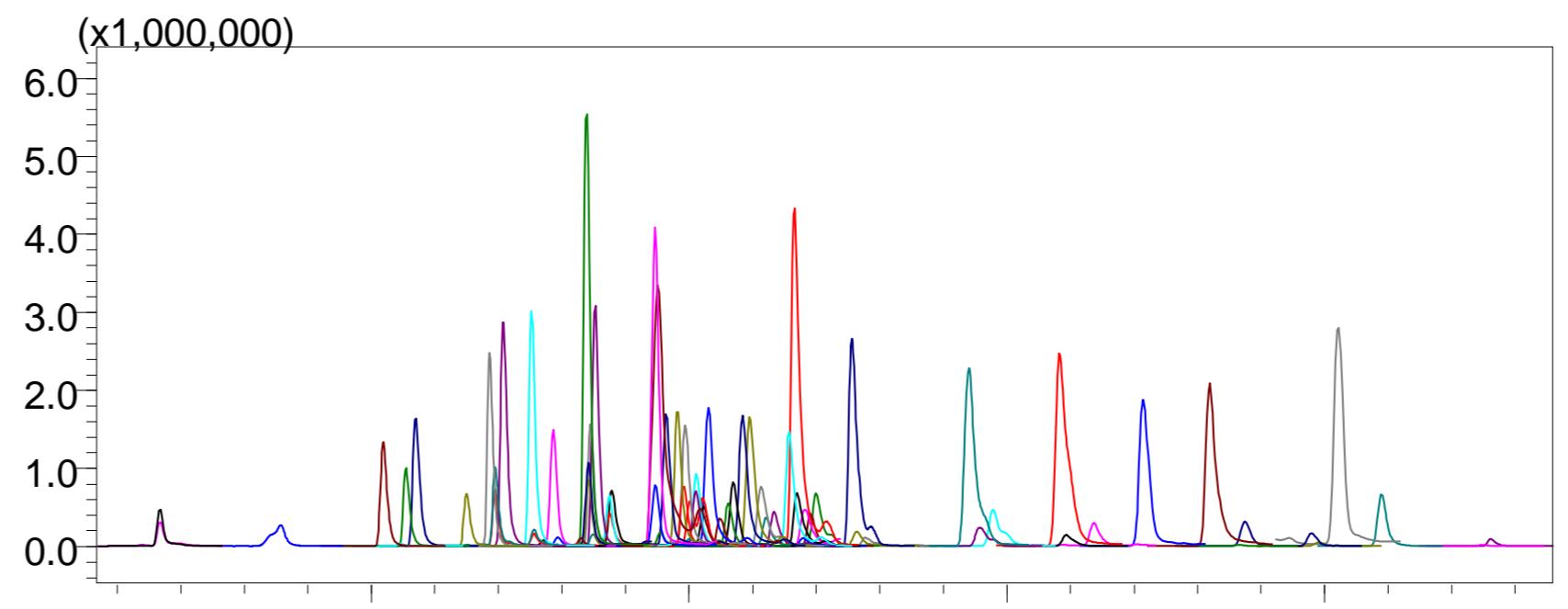


Figure 4 Representative Chromatogram for 84 pesticides using ESI-LCMS

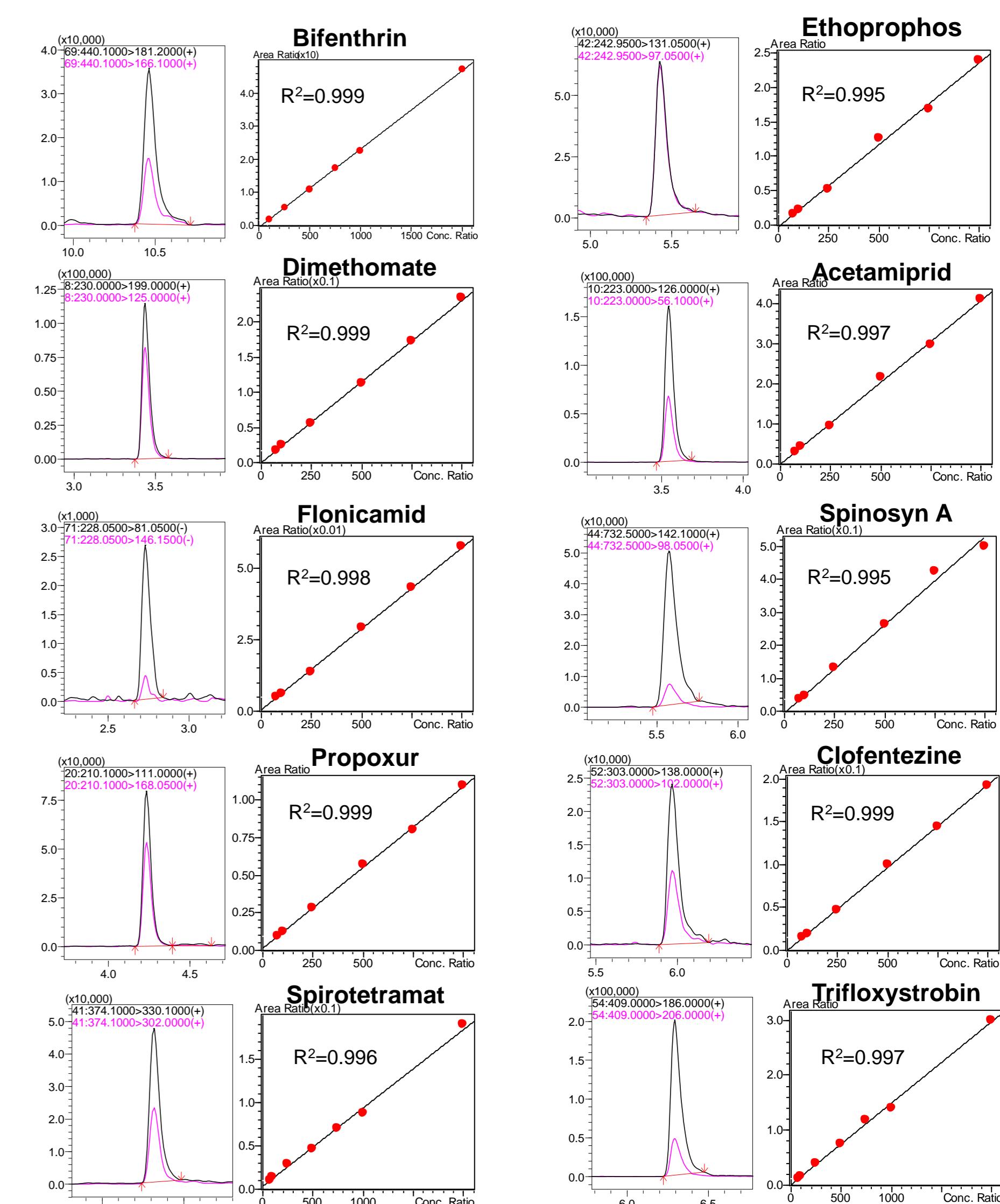


Figure 6 Calibration Curves and MS Chromatograms at 75 ng/g for ESI

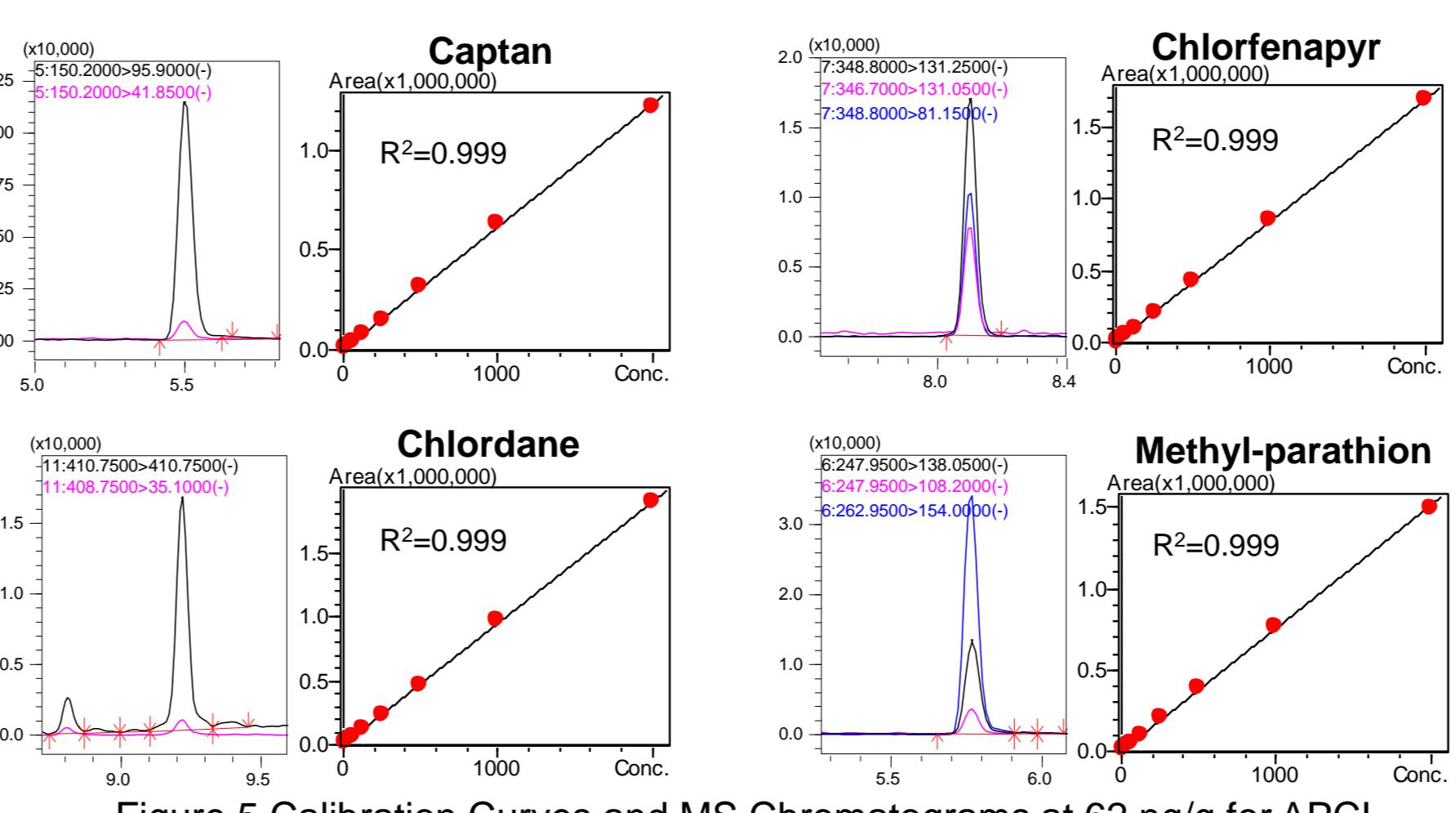


Figure 5 Calibration Curves and MS Chromatograms at 62 ng/g for APCI

4. Conclusions

A complete LCMS solution was developed for residual pesticide testing in cannabis matrix utilizing both APCI and ESI ionization techniques coupled with a single Shimadzu LCMS-8060. The APCI-LCMS method was developed and tested in cannabis flower matrix for the analysis of 10 California and Oregon regulated pesticides that have been traditionally analyzed by GCMS. An ESI-LCMS was further optimized and tested in cannabis flower matrix for the analysis of 84 total pesticides. The LOQs determined in this method were well below the action limits required by California and Oregon, demonstrating the viability of an LCMS total solution for cannabis testing in these two programs. The use of the ultrafast polarity switching capability of the LCMS-8060 allowed for accurate and sensitive quantitation of all 94 pesticides currently being regulated by California or Oregon.